

Data Sheet

HUMAN β -DEFENSIN 3 (aa 6-22)

ANTIBODY, MONOCLONAL

| | |
|----------------------------|---|
| Catalog no.: | AE1050.1 / AE1050.2 |
| Immunogen: | Synthetic human β -Defensin 3 (aa 6-22) (LQKYCRVRGGRCVLS) |
| Synonyms: | hBD-3, Beta-defensin 103, DEFB-3, Defensin-like protein |
| Swiss-Prot No: | P81534 |
| Gene Information: | Gene Name: DEFB103A, BD3, DEFB103, DEFB3 and DEFB103B GeneID: 55894 and 414325 |
| Host: | Mouse Balb/c |
| Clone no.: | L3-18b-E1 |
| Isotype: | IgG ₁ |
| Matrix: | Cell culture supernatant, Protein G purified, 50 mM TRIS pH 7.4 |
| Specificity: | Human β -Defensin 3 (aa 6-22), rat β -Defensin 3 There was no cross reactivity obtained with human β -Defensin 1, β -Defensin 2, and β -Defensin 4. |
| Contents: | 10 μ g / 100 μ g (lyophilized) Resuspend in 10 μ l / 100 μ l aqua bidest. |
| Known applications: | ELISA, RIA, Western Blot ¹ , immunohistochemistry (paraffin sections, 1:50) ¹ This antibody has not been tested for use in all applications. This does not necessarily exclude its use for non-tested procedures. The stated dilutions are recommendations only. We suggest that the applicant titrates the antibody in his/her system using appropriate negative/positive controls. |
| Store at: | 2-8 °C (lyophilized); - 20 °C (dissolved) Repeated thawing and freezing must be avoided |
| References: | 1. Wu Q, Gui P, Yao S, Zhu H, Li J, Li Y (2011) Expression of beta-Defensin-3 in Lungs of Immunocompetent Rats with Methicillin-Resistant Staphylococcus aureus Ventilator-Associated Pneumonia. <i>The Journal of surgical research</i> 169 (2): 277-283. |
| Last updated on: | 29 November 2024 |



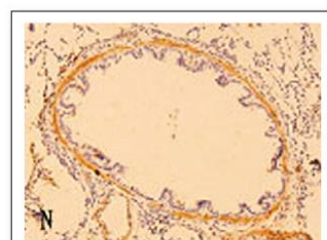


Figure 1: Immunohistochemistry image of β -Defensin 3 staining in paraffin section of rat lung tissue. Antigen retrieval was performed in 0.01 M sodium citrate buffer at 55–58°C for 5 min. The section was incubated with AE1050 (1:50) overnight and detected using a biotinylated secondary antibody in combination with a streptavidin-peroxidase conjugate. DAB was used as the chromogen. The section was counterstained with hematoxylin. AE1050 stains the bronchiolar epithelium in methicillin-resistant *Staphylococcus aureus*-inoculated rats. Original magnification: x200.

Wu Q et al. (2011) *J Surg Res* 169(2): 277-283

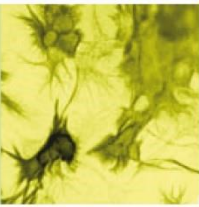

Antibodies

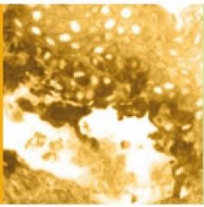

For research use only

Publishing research using AE1050? Please let us know so that we can cite your publication as a reference.



Immundiagnostik AG

Stubenwald-Allee 8a · 64625 Bensheim · Germany

Phone: +49 6251 70190-0 · Fax: +49 6251 70190-363 · dept.immunochemicals@immundiagnostik.com · www.immundiagnostik.com